

Jazzy® 1121

Tight-space indoor manoeuvrability is made easy with the Jazzy® 1121's compact overall size and mid-wheel drive design. Yet, the 1121 doesn't stop there. Outdoor performance—climbing and terrain transitions—is ensured through the use of Active-Trac® Suspension and the responsive action of the front anti-tips and rear casters. The Jazzy 1121 is the perfect power chair for your diverse and demanding lifestyle.



Jazzy® 1121
shown in Blue


21 Healey Rd., Dandenong, Victoria 3175
ACN 088 609 661 ABN 66 088 609 661
(Tel) 03 9706 4611 • (Fax) 03 9706 4622 • www.pridemobility.com

Pride®
Mobility Products
Australia Pty. Ltd.

SPECIFICATIONS

Weight capacity	136 kg (300 lbs.)
Maximum speed	Up to 8 km/h (5 mph)
Ground clearance	9 cm (3.5")
Turning radius	57 cm (22.5")
Overall length	98 cm (38.5")
Overall width	60 cm (23.8")
Wheels	Drive: 35.5 cm (14") aquatread
Front anti-tips	15 cm (6") solid
Rear castors	20 cm (8")
Suspension	Active-Trac®
Drive train	Two-motor, mid-wheel drive
Braking system	Intelligent braking (electronic, regenerative disc brakes)
Standard electronics	70 amp, PG VSI
Specialty controls	Yes
Battery charger	5 amp, on-board (standard)
Per-charge range	Up to 25 km (15.5 miles)
Battery requirements	(2) 12 volt, deep cycle
Battery size	55 AH (NF-22)
Battery weight	17 kg (38 lbs.) each
Base weight	61.5 kg (135.5 lbs.)
Standard seat weight	22.5 kg (49.6 lbs.) (euro-comfort)
Std. seat dimensions	43 cm x 43 cm (17" x 17") or 48 cm x 48 cm (19" x 19")
Warranty	2-year limited

COLOUR CHOICES

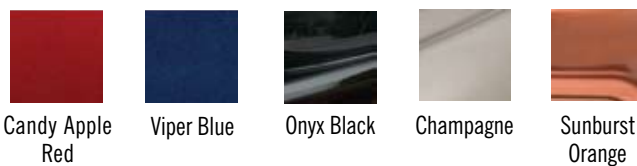

FEATURES

- Powerful, patented two-motor mid-wheel drive design provides superb power, acceleration and manoeuvrability
- Active-Trac® suspension for outstanding outdoor performance
- Climbs kerbs up to 7.5 cm (3")
- Adjustable front anti-tips to meet your preferences and terrain challenges
- Low centre of gravity design delivers a very smooth and stable ride
- Speed range adjustable up to 8 km/h (5 mph)
- State-of-the-art, programmable microprocessor-based joystick provides very smooth handling and control

OPTIONS

- Swing-away joystick
- Swing-away footrests
- Elevating leg rests
- Oxygen tank holder
- Cane/crutch holder
- Rear basket

Jazzy® 1121
shown in Blue



Pride®
Mobility Products
Australia Pty. Ltd.

21 Healey Rd., Dandenong, Victoria 3175

ACN 088 609 661 ABN 66 088 609 661

(Tel) 03 9706 4611 • (Fax) 03 9706 4622 • www.pridemobility.com

The information contained herein is correct at the time of publication; due to Pride's commitment to constant improvement and development, they reserve the right to alter specifications without prior notice. Theoretical range is calculated under test conditions in accordance with European standards. Actual range in normal use will depend on many factors, including the condition of the vehicle and its batteries, the weight of the driver, correct tyre pressure, ambient temperature, and the gradient and surface of the road or pavement.